

## Quick Reference Tools Guide

### Taking a Radiograph




**Add**  
Adds a new patient for radiographs



**Flip**  
Flip the radiograph vertically




**Save + Acquire**  
Saves current patient information and begins the acquisition process

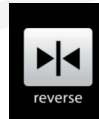

**Reverse**  
Flip the radiograph horizontally




**Anatomy**  
Choose a different anatomy /view

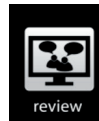

**Review**  
Review the radiographs taken in the current study



**Views**  
Choose a different view

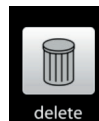

**Delete**  
Delete a study of patient



**Left**  
Rotate the radiograph 90° left



**Save**  
Saves any positioning changes made to the selected radiograph (left, right, flip, reverse)



**Right**  
Rotate the radiograph 90° right




**Patient**  
Returns to the Patient Screen. If done from the Review Screen, it will send your x-rays to a DICOM server

### Adding Images to an Existing Study



**Acquire**  
Add images

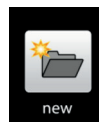

**New**  
Create a new study



YOUR IMAGE IS OUR BUSINESS.



a VCA ANTECH company

800.268.5354 [soundeklin.com](http://soundeklin.com)